

# ***EC-720R/M,I,L***

***Lower G.I. tract scopes***




# EC-720R/M,I,L

## Lower G.I. tract scopes



### Wide 170° field of view

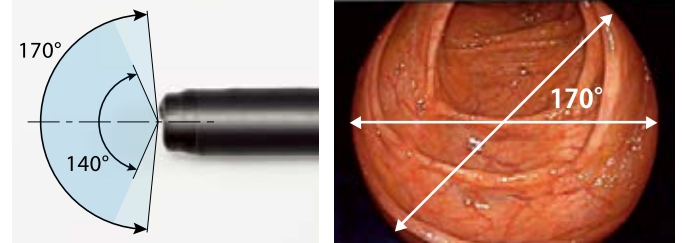

Even areas that are hard to observe such as the reverse side of folds could be observed and approached smoothly.

### Advanced Force Transmission

The flexible portion is designed to transmit operator's movements, pushing and rotating, to the distal end of endoscope.




### Adaptive Bending

The end of flexible portion is soft, allowing the scope to bend easily. Flexible portion is elastic, and easy to return to its straight shape.






### Insertion Performance

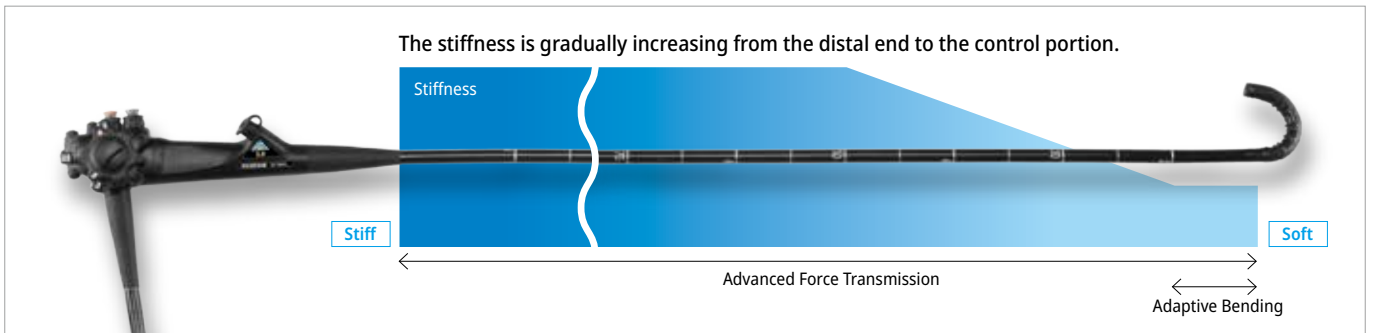

### Specification

EC-720R/M,I,L	
Field of view	170°
Viewing direction	0° (Forward)
Observation range	2~100 mm
Bending capability	UP: 180° DOWN: 180° RIGHT: 160° LEFT: 160°
Working length	1,330 mm(M)/1,520 mm(I)/1,690 mm(L)
Total length	1,630 mm(M)/1,820 mm(I)/1,990 mm(L)
Distal end diameter	12.8 mm
Insertion tube diameter	12.8 mm
Minimum instrument channel diameter	3.8 mm
Image size	Super image

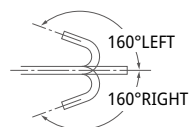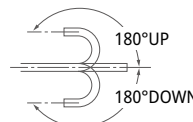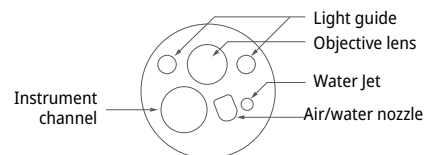

Product name: Video Endoscope GMDN: 36117  
Generic name: Flexible video colonoscope, reusable