


HANAR STENT[®]

Esophagus Skidproof (CCC)

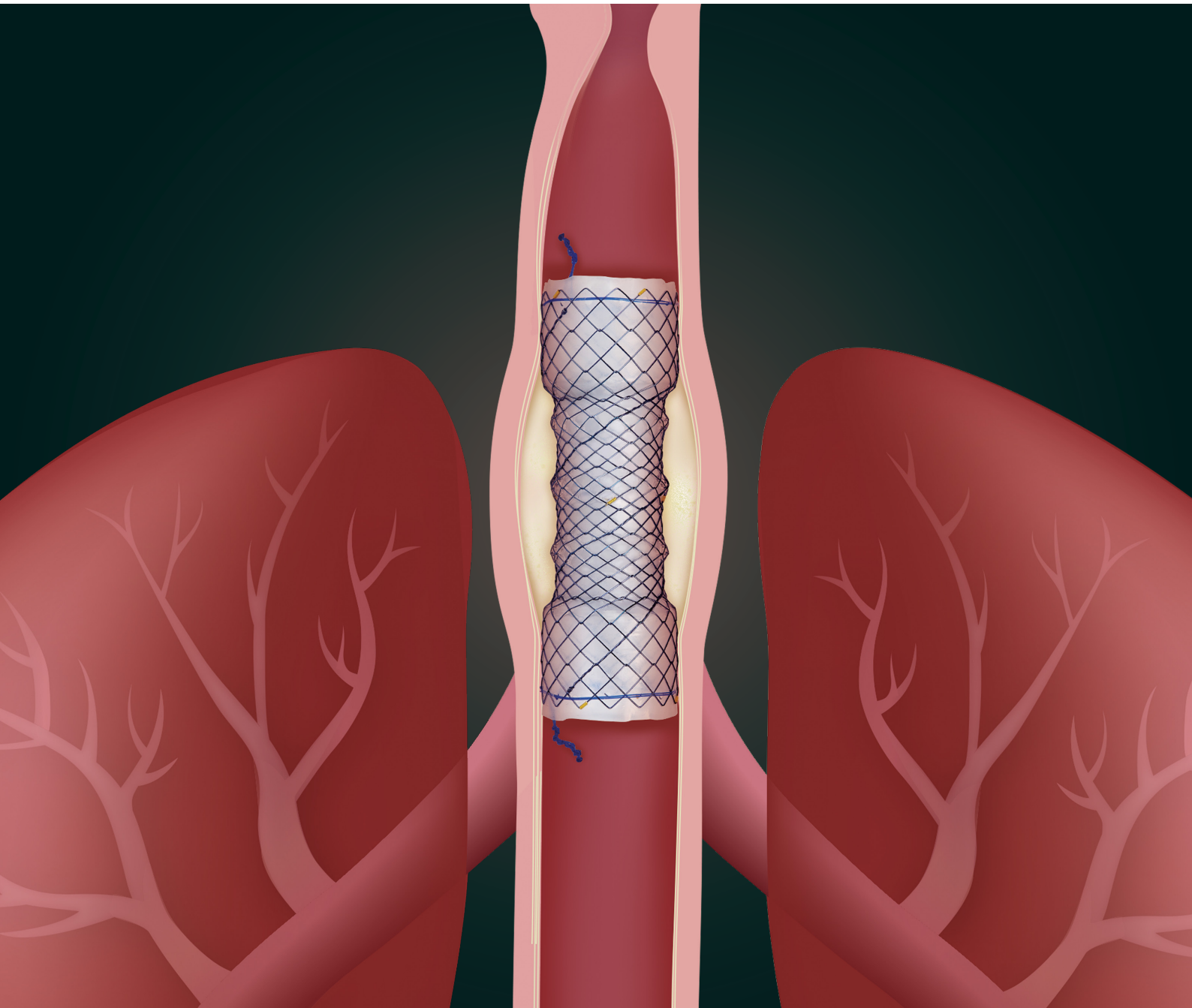


Skidproof body structure for Anti-Migration.

- Fully covered stent


HANAR STENT®

Esophagus Skidproof (CCC)

Features

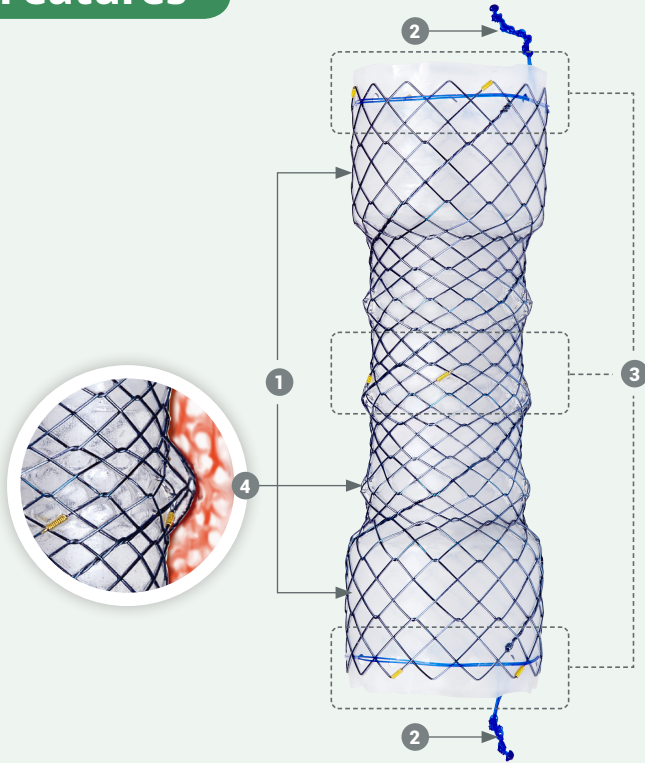

- 1 **Dumbbell heads on the both ends for Anti-Migration**
- 2 **Lassos on the both ends for repositioning**
- 3 **Radiopaque markers for precise stent placement**
- 4 **Skidproof body structure for Anti-Migration**
Multiple protuberances on its body helps to resist migration.
- 5 **Hook&Cross Weaving Structure**

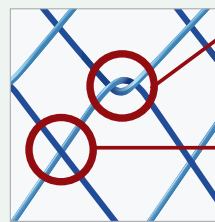

Hook

Flexibility to conform within the anatomy

Cross

High pushability through catheter during deployment within the anatomy

Ordering Information

※ Sizing and availability varies by country

● ECD

Model	Stent (mm)			Delivery Device	
	Diameter	Usable Length*	Total Length*	Length (mm)	Diameter (mm/Fr)
ECD-20-080-070	26-20-26	40	80	700	6/18
ECD-20-100-070		60	100	700	6/18
ECD-20-120-070		80	120	700	6/18
ECD-20-140-070		100	140	700	6/18
ECD-20-160-070		120	160	700	6/18

* Increment of 20mm